


Mini Pillow



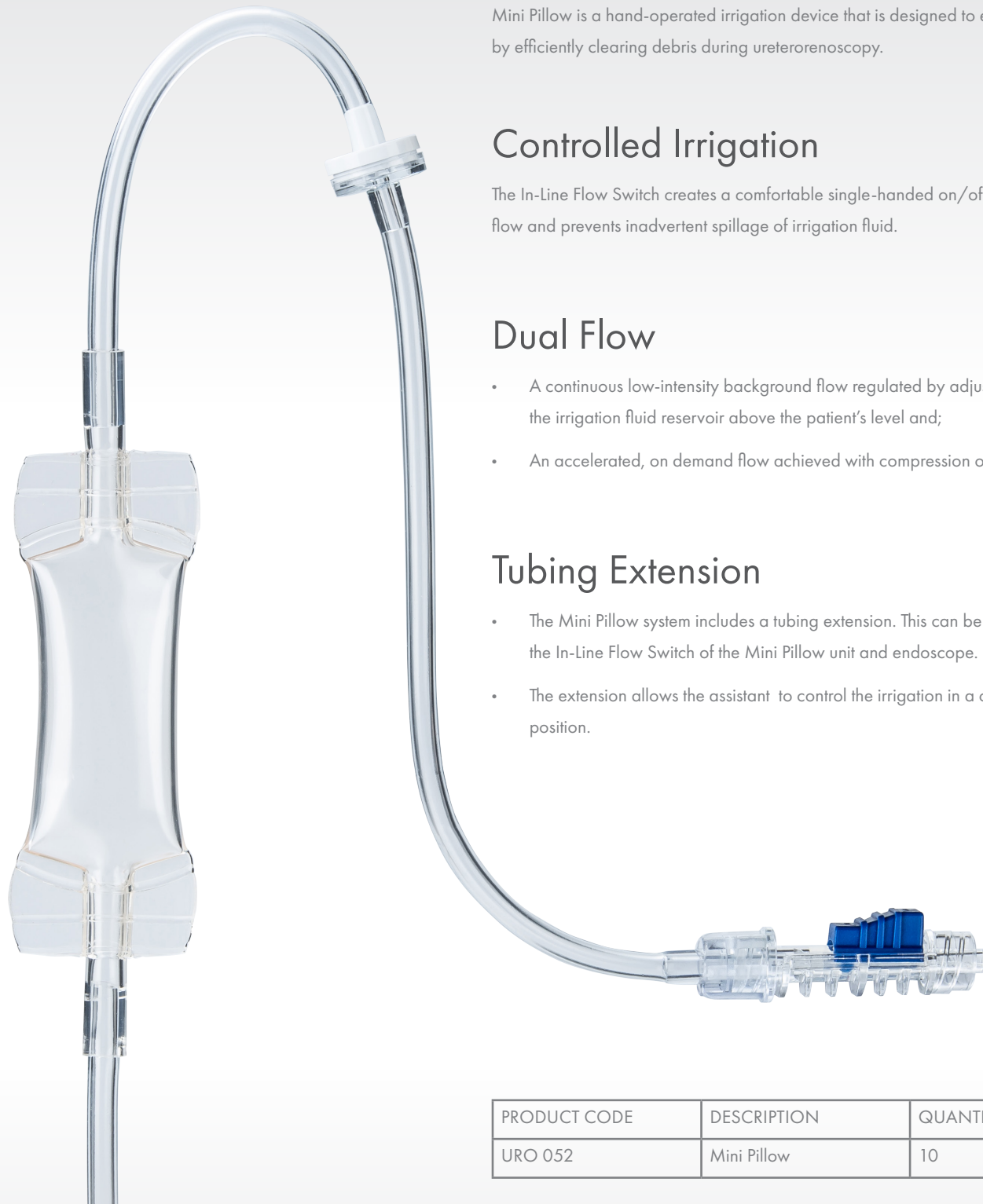

Mini Pillow

Mini Pillow is a hand-operated irrigation device that is designed to enhance visibility by efficiently clearing debris during ureterorenoscopy.

Controlled Irrigation

The In-Line Flow Switch creates a comfortable single-handed on/off control of the flow and prevents inadvertent spillage of irrigation fluid.

Dual Flow

- A continuous low-intensity background flow regulated by adjusting the height of the irrigation fluid reservoir above the patient's level and;
- An accelerated, on demand flow achieved with compression of the pillow.

Tubing Extension

- The Mini Pillow system includes a tubing extension. This can be attached between the In-Line Flow Switch of the Mini Pillow unit and endoscope.
- The extension allows the assistant to control the irrigation in a comfortable position.

PRODUCT CODE	DESCRIPTION	QUANTITY
URO 052	Mini Pillow	10



Wismed CC
71 Princess Alice Avenue,
Durban, 4001, South Africa
Tel: +27 31 206 0177
Email: orders@wismed.co.za
www.wismed.co.za



TNMC Devices Limited (Llc)
Krežmina 1
10000 Zagreb, Croatia
Tel: +38 59 1622 9658
andretn@tnmc.biz

